



# T.E.D.\* Anti-Embolism Stockings

## Latex-Free

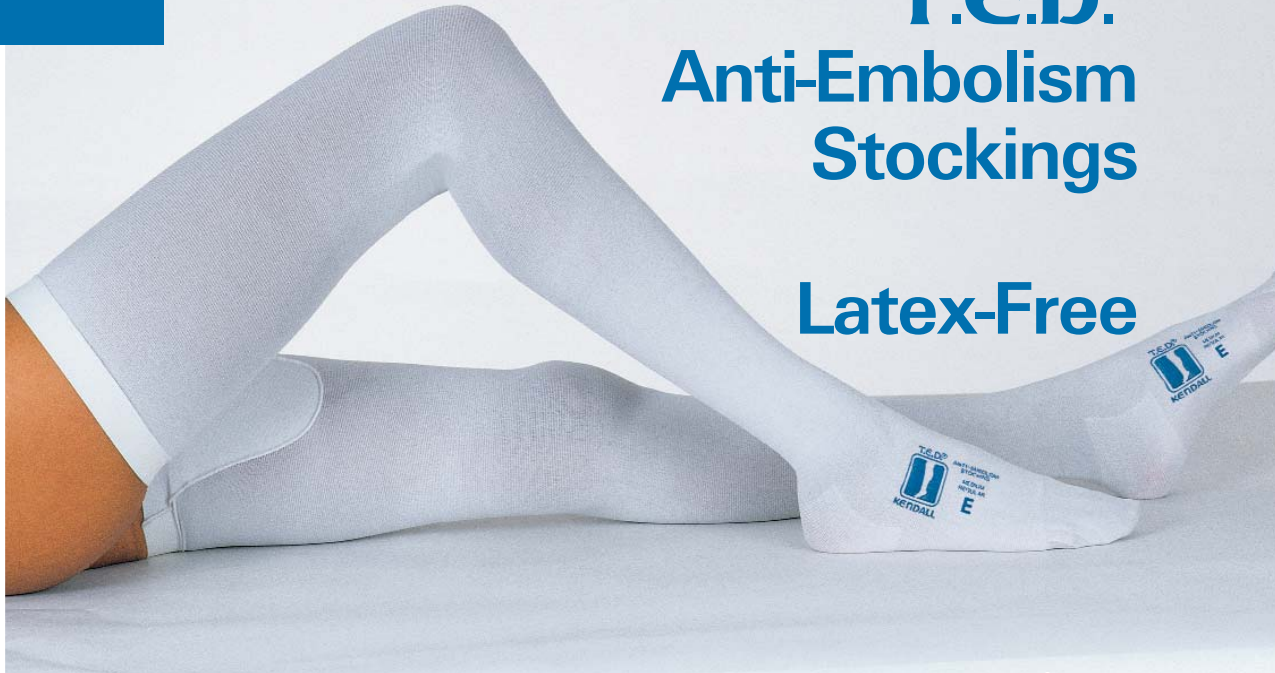

If proven efficacy in an anti-embolism stocking is important to you, be aware that not all stockings are the same. Only the T.E.D.\* Anti-Embolism Stocking has been clinically proven to prevent DVT in over 14,000 patients.

## The new T.E.D.\* unique silicone coated thigh band

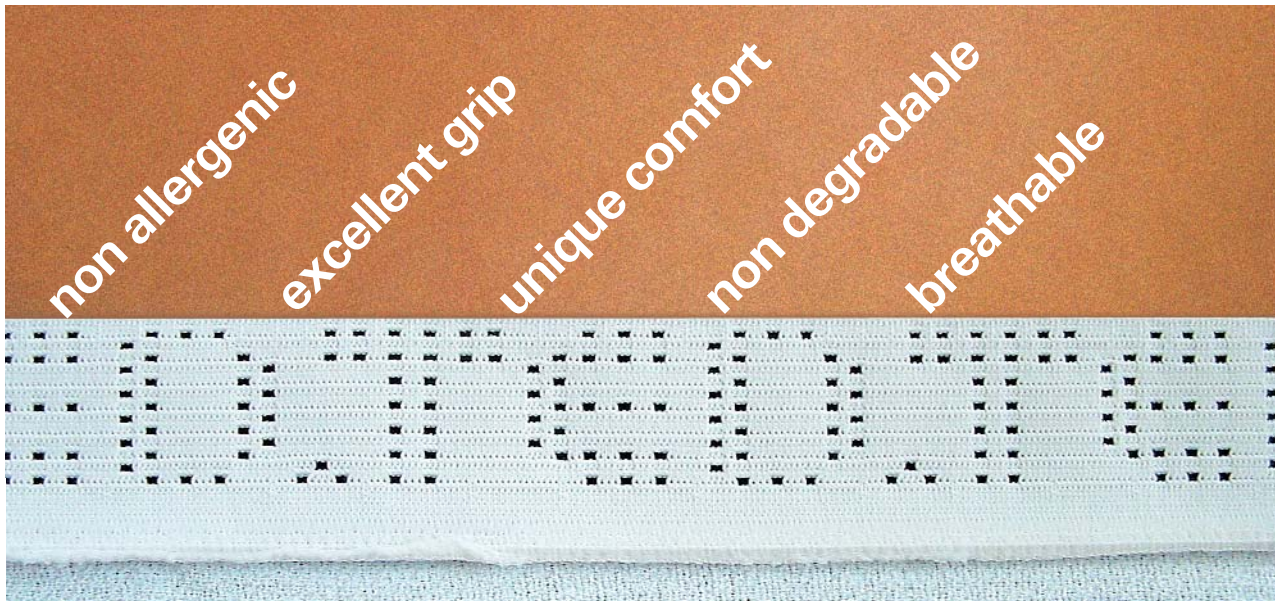

Patient Comfort  
Assured Efficacy  
Safety by Design  
Ease of Use

To further improve our highest quality **T.E.D.\* Anti-Embolism Stockings**, we have developed a flexible film that is wider and thinner than most other available coatings, spreading the forces over a greater area and therefore

maximising the performance and comfort of the product. When a product has to stay in place like the **T.E.D.\* Anti-Embolism Stockings**, the use of a flexible thin silicone elastomeric is an ideal way to provide the

non-slip properties that are necessary. The **T.E.D.\*** silicone coating is comfortable against the skin whilst ensuring the stockings stay in place.

**tyco**

Healthcare

**KENDALL**

Latex-Free  
NEW!

# T.E.D.\* Anti-Embolism Stockings

## All Styles 100% Latex Free

**Count on T.E.D.\* Stockings to:  
Reduce DVT by 50% – Prevent venous distension – Prevent venous stasis**

### Assured Efficacy

- **Reduce DVT by 50%<sup>(1)</sup>**  
Proven in a wide variety of medical/surgical specialities.
- **Improve patient outcomes**  
Combining the T.E.D.\* Stocking with mechanical or pharmaceutical prophylaxis will lead to a lower incidence of DVT, better patient outcomes and lower healthcare costs.
- **Provide baseline DVT protection**  
Whatever the patient's level of risk or when pneumatic compression is interrupted, T.E.D.\* Stockings will help to ensure uninterrupted baseline protection.
- **Protect after discharge**  
T.E.D.\* Stockings worn post-discharge can help reduce the chance of delayed DVT complications.

<sup>(1)</sup> References available upon request.

### Safety by Design

- **100% latex free available for all styles**
- **New silicone band**  
for excellent grip and unique patient comfort due to non-allergic material and increased breathability.
- **A clinically proven pressure pattern**  
of 18 mmHg at the ankle, 14mmHg at the calf, 8 mmHg at the knee, 10 mmHg at the mid-thigh, 8 mmHg at the upper thigh increases blood flow velocity to 138% and reduces venous distension.
- **Inlay circumferential knit**  
with one-way stretch delivers the clinically proven pressure pattern to the appropriate areas of the leg.
- **A pressure break at the popliteal vein**  
ensures that blood will continue to flow smoothly through this critical area.
- **Interrupted band and 2-ply panel**  
prevent a tourniquet effect at the femoral vein.

### Ease of Use

- **Fits up to 95% of all patients**  
T.E.D.\* Stockings come in three styles – Knee Length, Thigh Length and Thigh Length with Belt – with 27 different sizes to accommodate most patients.
- **A defined heel pocket**  
ensures correct placement.
- **Inspection toe opening**  
allows for quick inspection of the skin and pulse.
- **Sophisticated manufacturing techniques**  
and strict quality control help to ensure that the correct pressures will be delivered and the stockings will be free of defects.

### T.E.D.\* thigh length stockings



REF	calf size	length
3071LF	small	short
3130LF	small	regular
3222LF	small	long
3310LF	medium	short
3416LF	medium	regular
3549LF	medium	long
3634LF	large	short
3728LF	large	regular
3856LF	large	long

### T.E.D.\* thigh length with belt



REF	calf size	length
3306	extra-small	regular
3320	extra-small	long
3039	small	regular
3364	small	long
3144	medium	regular
3449	medium	long
3221	large	regular
3523	large	long
3922	extra-large	regular
3995	extra-large	long

### T.E.D.\* knee length stockings



REF	calf size	length
7071	small	regular
7339	small	long
7115	medium	regular
7480	medium	long
7203	large	regular
7594	large	long
7604	extra-large	regular
7802	extra-large	long

\*Trademark © 2003 Tyco Healthcare Group LP. All rights reserved.

M-TED LF/ GB

**tyco**  
Healthcare

**KENDALL**

### For more information please contact:

**Tyco Healthcare UK Limited**  
154 Fareham Road  
GOSPORT  
Hampshire  
PO13 0AS  
UK  
Tel: +44 (0) 1329 224226  
Fax: +44 (0) 1329 224334

**Tyco Healthcare**  
Birchwood Court  
Montrose Road - P.O. Box 8108  
Halfway House  
1685 MIDRAND  
SOUTH AFRICA  
Tel: +27 (0) 11 265 2000  
Fax: +27 (0) 11 265 2050

**Tyco Healthcare Regional Export**  
154 Fareham Road  
GOSPORT  
Hampshire  
PO13 0AS  
UK  
Tel: +44 (0) 1329 224258  
Fax: +44 (0) 1329 224086

**Tyco Healthcare Norden AB**  
Gårdsvägen 18  
SE-169 67 SOLNA  
SWEDEN  
Tel: +46 (0) 8 585 605 05  
Fax: +46 (0) 8 585 605 06  
www.tycohealth.com

**Tyco Healthcare Nederland BV**  
Hogeweg 105  
5301 LL ZALTBOMMEL  
The NETHERLANDS  
Tel: +31 418 57 66 00  
Fax: +31 418 57 67 93  
TycoHealthNed@TycoHealth.com  
www.tycohealthcare.nl

**Tyco Healthcare Belgium N.V.-S.A.**  
Generaal de Wittelaan 9/5  
B-2800 MECHELEN  
BELGIUM  
Tel: +32 (0) 15 29 44 50  
Fax: +32 (0) 15 29 44 55  
tycohealthcare.belgium@tycohealth.be

04/2003